
## **Packing Specifications**

Please retain the original packaging and use it whenever storing or transporting the scale. Always pack the scale as shown below. When packing the scale, be sure to remove the pan and the pan base. Shipping the scale with the pan installed may damage the unit. Always place the items in their designated compartments, as shown below.

## Scale



Do not pack the scanner or any scanner components with the scale. Always use a separate box designated specifically for the scanner. Before packing the scanner disassemble the stand by separating the Apron from the stand. To remove the Apron, unscrew the two M3 x 6mm screws as shown in the picture below.

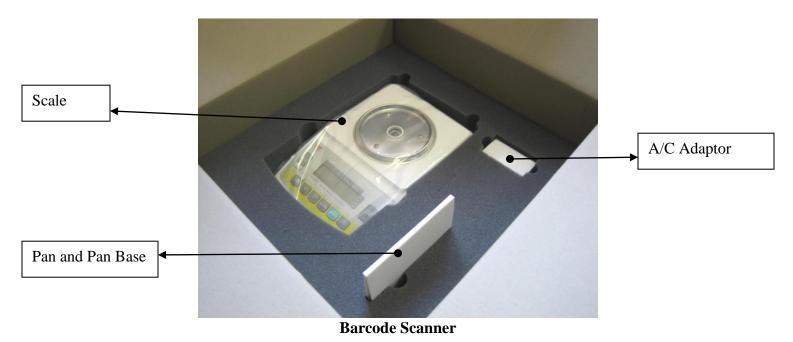